

TSIII ø5.5 Implant

OSSTEM[®]
IMPLANT


Specification of ø5.5 Fixture




System	TSIII
Diameter	Ø 5.5
Length	6, 7, 8.5, 10, 11, 13, 15mm
Surface	SA / CA
Mount	Pre / No

Indication of ø5.5 Fixture


1. To overcome Poor Initial Fixation of ø5.0 Fixture : Replace 5.0 Fixture with 5.5 Fixture when the initial fixation is weak

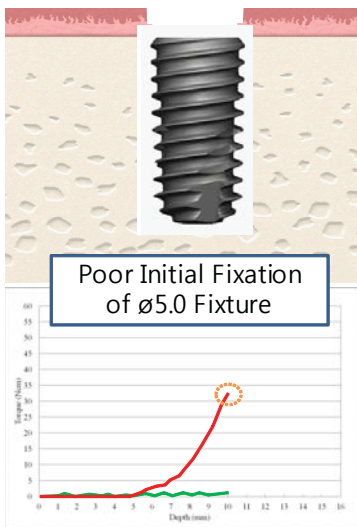

2. Immediate Replacement for Failed ø4.5 Fixture : When ø4.5 Fixture failed and bone quality is weak, ø5.5 Fixture is a good alternative for the position

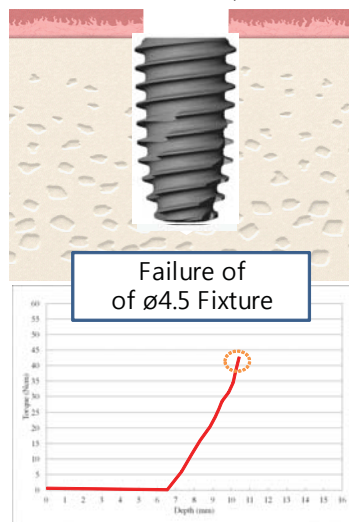

3. Not enough Ridge for Ultra-Wide Fixture : When there is adjacent teeth or not enough ridge for Ultra-Wide fixture, ø5.5 will be selected rather than ø6.0

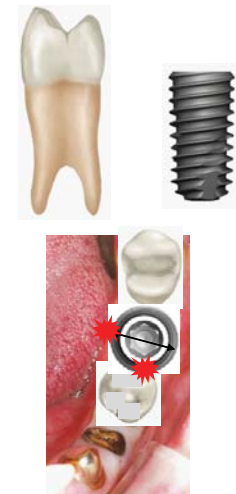


OSSTEM Fixture Line-up

	Ø3.0	Ø3.5	Ø4.0	Ø4.5	Ø5.0	Ø5.5 New	Ø6.0	Ø7.0
TS								
US	-							
SS	-							

Ø5.5 * 10mm




Hard
Bone

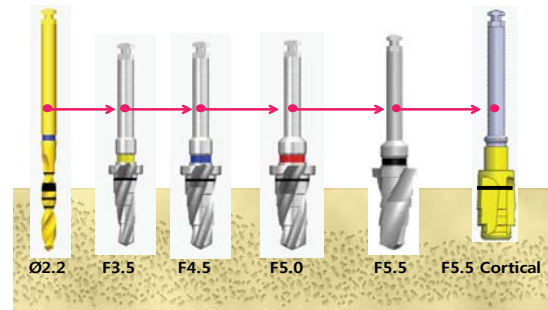


○ → ○ → ○ → ○ → ○ → ○ Drilling 6 times

Normal
Bone

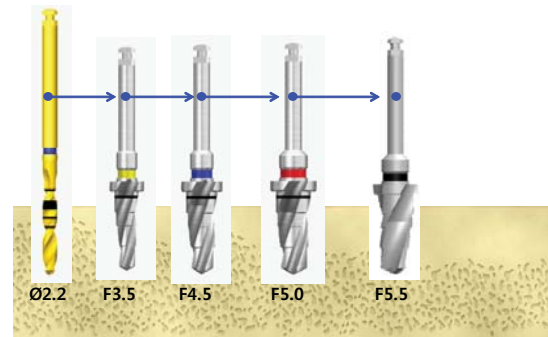

○ → ○ → ○ → ○ → ○ Drilling 5 times

Soft
Bone

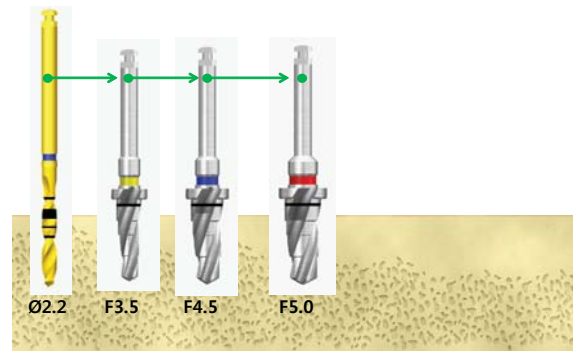

○ → ○ → ○ → ○ Drilling 4 times

OSSTEM[®]
IMPLANT

DEUTSCHE OSSTEM GmbH

Mergenthalerallee 25, 65760 Eschborn, Germany

Tel: +49(0)6196 777 550 Fax: +49(0)6196 777 5529

E-mail: info@osstem.de www.osstem.de